



Cattani air care


AC400

CATTANI
air technology

OIL-LESS COMPRESSOR

Il compressore a secco Cattani (privo di lubrificazione ad olio) completo di impianto di essiccazione e di una serie di filtri, è in grado di fornire aria compressa secca e di elevata qualità e purezza. Con il filtro sterilizzante (opzionale) l'aria viene inoltre trattata con efficienza del 99,9999+% su particelle da 0,01 µm. Abbiamo fatto della qualità ed affidabilità la nostra bandiera, per questo siamo in grado di dare 3 anni di garanzia.

The CATTANI oil-less compressor (without oil lubrication) complete with drying system and a set of filters, is capable of supplying dry compressed air of extremely high quality and purity. With the sterilising filter air (optional) is treated with an efficiency of 99.9999+% on 0.01µm particles. We have chosen quality and reliability as our emblems: this is why we can grant you a 3-year warranty.

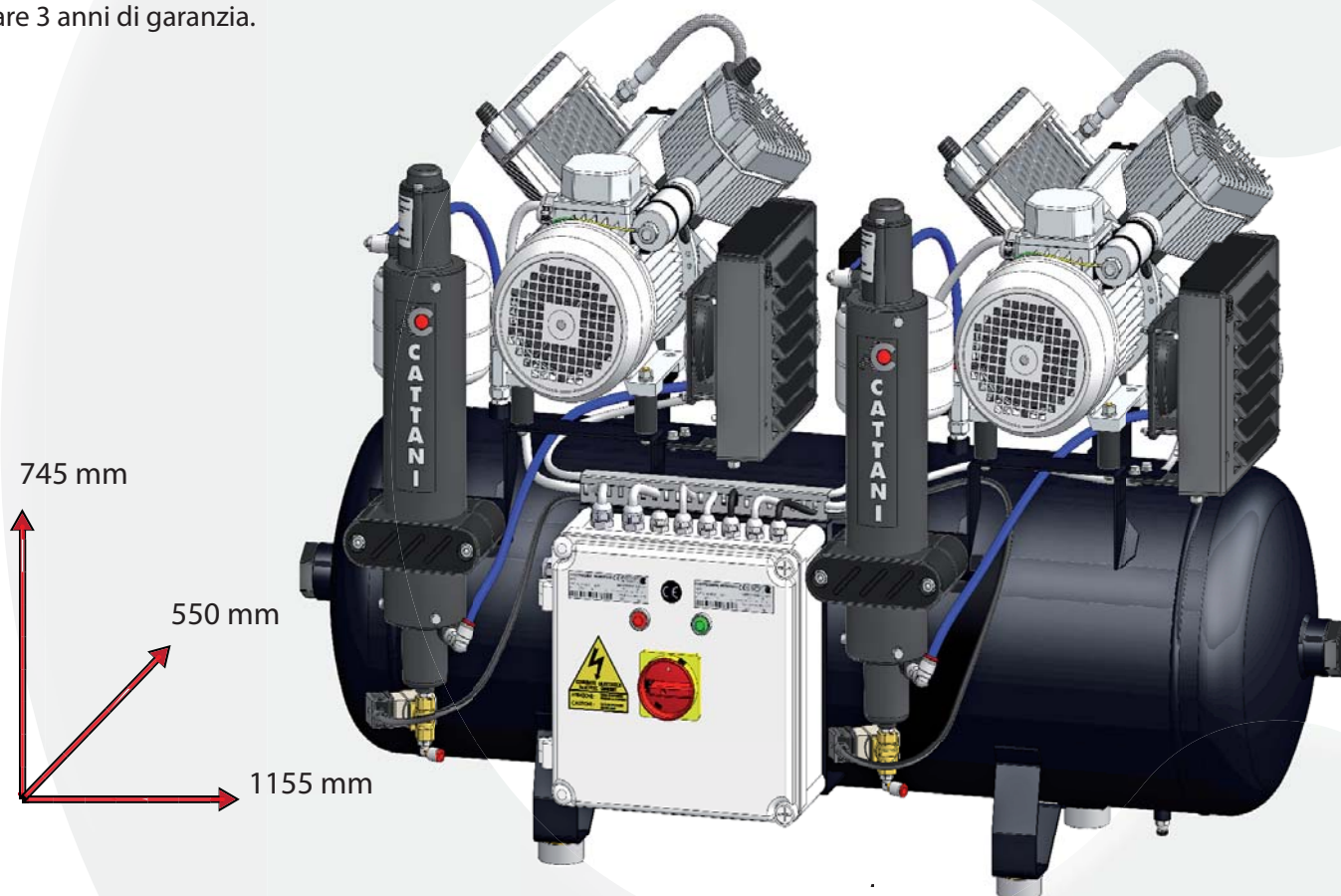

Tipo	Dati elettrici	Codice
AC400	(1~) 230-50 hz- 15,4 A - 2,4 kw	013430
AC400	(3~) 400-50 hz- 7,2 A - 3 kw	013435
AC400	(1~) 220-60 hz- 18,4 A - 3 kw	013407

Aria prodotta NI/min at 5 bar	320
Capacità serbatoio litri	100
Peso kg	112
Livello pressione sonora db (A)	73

Type	Electrical data	Code
AC400	(1~) 230-50 hz- 15,4 A - 2,4 kw	013430
AC400	(3~) 400-50 hz- 7,2 A - 3 kw	013435
AC400	(1~) 220-60 hz- 18,4 A - 3 kw	013407

Delivered air NI/min at 5 bar	320
Tank capacity liters	100
Weight kg	112
Sound pressure level db (A)	73

CATTANI S.p.A.

Via Natta, 6/A – 43122 Parma – Italy – Tel. +39 0521 607604
 SALE DEPT. FAX: +39 0521 607628 – PURCHASING DEPT. FAX: +39 0521 607855
 ACCOUNTING DEPT. FAX : +39 0521 399966
<http://www.cattani.it> Email: info@cattani.it PEC: cattani@peclife.it

Codice Fiscale e Partita I.V.A. 01720020344 – E.E.C. VAT IT 01720020344 –
 Capitale Sociale € 1.549.800,00 I.V. - R.E.A. 173616
 Registro Imprese Parma n. 01720020344

AZIENDA CON SISTEMA DI GESTIONE
 PER LA QUALITÀ CERTIFICATO DA DNV
 = UNI EN ISO 9001:2008 =